

HOW TO USE YOUR BREZTRITM AEROSPHERE[®] INHALER

Important Information

- For oral inhalation use only.
- Use BREZTRITM AEROSPHERE[®] exactly as your doctor tells you to.
- If you have any questions about the use of your inhaler, ask your doctor or pharmacist.
- Clean your inhaler once each week. It is very important to keep the plastic actuator clean so that medicine will not build up and block the spray through the mouthpiece.

Your BREZTRITM AEROSPHERE[®] inhaler should be stored at room temperature.


Priming your inhaler

Get started with these steps:


1. Remove cap, and check inside the mouthpiece for objects.




3. Press the centre of the indicator (which tracks dosing) until it cannot be pressed further—a puff of medicine will be released from the mouthpiece (you may hear a soft click). Repeat three more times (4 times total), shaking before each puff.

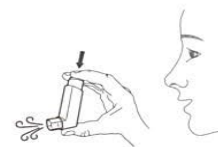

2. Hold your inhaler upright, away from your face. Shake well.

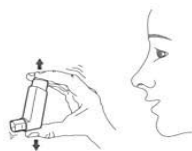

4. The dose indicator should now point to the right of "120". This means your inhaler is ready and primed.

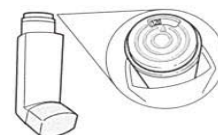

HOW TO DO WEEKLY RINSE FOR YOUR BRETZI AEROSPHERE INHALER

Weekly rinse


Rinse the **yellow actuator** weekly so that the medicine does not build up and block the spray through the mouthpiece. It's a good idea to mark this day in your calendar or add a reminder on your phone.

- Rinse the yellow actuator weekly.
- Do not allow the canister to get wet.
- Re-prime after rinsing.

1. Remove the canister and set it aside. **Do not allow the canister to get wet.**



5. Look into the actuator and mouthpiece for medicine build-up. If there is any build-up, repeat Steps 3 through 5.

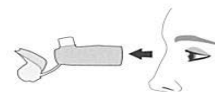

2. Remove the mouthpiece cover.



6. Air-dry, preferably overnight. **Do not** put the canister back into the actuator if it is still wet.



3. Run warm water through the mouthpiece for 30 seconds and then through the top of the actuator for 30 seconds. **Rinse for 60 seconds in total.**




7. When dry, **replace the mouthpiece cover first** and then gently press the canister down into the actuator.



4. Shake off as much water as you can. **Do not** dry with a towel or tissue.



8. Re-prime the inhaler by spraying **2 test-puffs**, shaking before each test-puff.

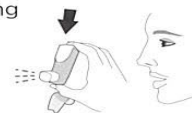

x2 total shakes and test-puffs

