

DRY POWDER INHALER

BREEZHALER



1

PULL OFF CAP



2

FLIP MOUTHPIECE TO OPEN



3

INSERT CAPSULE



4

FLIP MOUTHPIECE TO CLOSE



5

PRESS AND LET GO



6

BREATHE OUT




7

FAST DEEP BREATH IN & HOLD



8

REMINDER

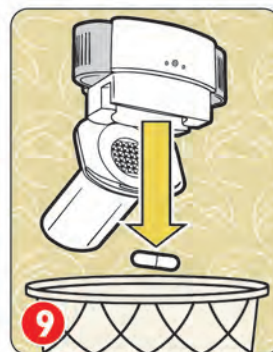

9

THROW AWAY

IMPORTANT INFORMATION

Expired?: Check the expiry date on the box or on the pack that holds the capsules.

Reminders: Keep the Breezhaler dry. Clean with dry cloth only.

Do not breathe into the device.

Capsules are sensitive to light, moisture and extreme temperatures.

Make sure to remove only **ONE** capsule from the pack at a time.

Do not touch **USED** capsules.

If the dry powder gets on your hands, make sure you wash your hands thoroughly.

Breezhaler device is recyclable.

This material is designed for information purposes only. It should not be used in place of medical advice, instruction and/or treatment. If you have specific questions, please consult your doctor or appropriate healthcare professional.

Illustrations developed by the Community Pediatric Asthma Service, www.ucalgary.ca/icancontrolasthma

For more information on inhaled medications and device demos, visit our website at:
www.ucalgary.ca/asthma

